FT 4925

Bedienungsanleitung/ Garantie

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Mode d'emploi

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Instruction Manual

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Product Overview

Thank you very much for choosing this product. This thermometer takes temperature in seconds by measuring the heat generated by the ear or the surface skin of forehead. It has the following advantages:

- 1) 4-in-1 design
Measures the ambient temperature and the body temperature in the ear or on the forehead. Furthermore, the current time is displayed.
- 2) Instant measurement
It provides you the reading in seconds.
- 3) Convenient, probe cover free and easy to use.
- 4) High speed accurate scan provides you precise reading of temperature.
- 5) Automatically converts forehead temperature display into ear temperature display with the removal of the front cover. Simple to use.
- 6) Alarm of fever
When the temperature reading exceeds 38.0°C, the thermometer will send out "Beep-Beep-Beep" for a long one followed by two short ones, and the "Fever" icon will go on .
- 7) Automatic display of memory
Every time you turn it on, the previous reading will be shown.

Icons on the Unit


Read the instructions carefully before use and keep them for your future reference!

Medical device type B

Precautions for Using the Product

- This thermometer is intended for the practical use at home.
⚠ It does not replace consulting a physician!
- With a waterproof tip, after each use, the thermometer can be cleaned with clean cloth or paper or alcohol. But the whole thermometer should not be immersed into any liquid.
- Do not expose the device to the extreme temperatures for a long time. For instance, in the car, under direct sunlight.
- Do not drop it or strike it with something else. Do not touch the lens of the probe with your fingers. Do not disassemble the device by yourself.
- The device is designed without ear caps , which are not necessary. Keep the probe clean, and protect it from earwax, sweat and oil stains. Otherwise, the infrared measurement function will be affected.
- Do not measure body temperature within 30 minutes after meals, exercise or bath.

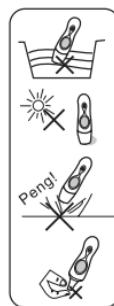- When the environment changes a lot (for instance, from indoor to outdoor), please wait in the new environment for about 30 minutes before measuring the new temperature.
- The intervals shall be 10 seconds when constant measurement are made. It is normal that there are deviations in constant measurement between short intervals, in which case the average reading can be taken as the result.
- There is no absolutely standard value of body temperature. For healthy persons the temperature can fluctuate between $0.2 \sim 1^\circ\text{C}$ ($32.36 \sim 33.8^\circ\text{F}$) in different areas of the body. Always measure the body temperature at the same time of the day.
- When the environmental temperature exceeds the functional temperature range of the device, the device will show "Err" symbol.
- When the temperature of the measured is lower than 32°C , "Lo" symbol will be displayed on the screen; when the temperature of the measured is higher than 42.9°C , "Hi" symbol will be displayed on the screen.

Location of Controls

- 1 Power / Scan button
- 2 Display
- 3 Front cover / Attachment for forehead measurement
- 4 Stand
- 5 Probe
- 6 Battery compartment
- 7 Time set button

Definition of Symbols

Reading Display

Measuring in Progress

Ear Scan Mode

Forehead Scan Mode

Fever

Celsius Scale

Fahrenheit Scale

Clock

Low Battery

Last Memorized Temperature Display

Description of LCD Display

Display	Meaning	Action
	Display the time and environment temperature in turns while not used to measure the baby temperature.	Set the time when you use the device for the first time or when you re-install the battery.
	Is displayed for 1 second when the Thermometer is first turned on.	Wait for the last memoried temperature to appear.
	Is displayed after the full display appears. Identifies the last measurement stored in the memory.	The thermometer is ready for measurement.
	The quantity of electrical charge is not sufficient.	Though the battery can be still used, you'd better change it as soon as possible.
	The battery is running low.	Unless you change the battery, you can't use the product to conduct measurement.
	Measured temperature too high.	When the temperature of the measured object is higher than 42.9°C.
	Measured temperature too low.	When the temperature of the measured object is lower than 32°C.
	Ambient temperature too high or too low.	When the environment temperature exceeds the working temperature scope of the device.

Display	Meaning	Action
	Measurement completed	The thermometer switches off automatically.

Time and Ambient Temperature

The device is supplied with an inserted battery. Remove the protective film between the battery and the terminal. This protects the battery from discharging during transport.

Set the time when you use the device for the first time or when you re-install the battery. Then, the function of displaying the environmental temperature and time will automatically start up.


1) Method of setting

Press the power button to switch the thermometer on.

Open the battery cap on the back, press the time-setting button to start hour (hour number flashes) and minute (minute number flashes) setting, and press the button of scan to adjust the numbers. Press the time set button to confirm (Refer to the Picture 6.1, 6.2).

2) Mode of display

When time has been set, the thermometer will display the time and environmental temperature in turns while not used to measure body temperature (Refer to the Picture 6.3, 6.4).

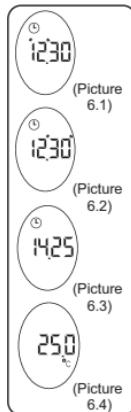Changing the Measurement Scale

The infrared thermometer can display temperature measurements in either °C or °F.

The operation how to change the measurement scale between °C and °F as: After all icons displaying and when the picture "--°C" or "°F--" displaying, please press the scan button, turn the °C or °F to °F or °C. Wait for 3 seconds and the thermometer will automatically enter the measurement mode.

Note:

Your thermometer provides both Fahrenheit and Celsius readings in countries where the Fahrenheit units are adopted, otherwise it provides only the Celsius unit. It isn't a bad one if your product can't change °C to °F.

Taking Temperature in Ear Temperature Mode

- 1) When the system is in the state of time/environmental temperature display, press the button of scan, and the last temperature measurement will be displayed; then after "Beep-Beep", the system enters into the mode of measurement.
- 2) Remove the probe cap, and the thermometer enters into the ear measuring mode (Refer to the Picture 7.1).
- 3) Carefully pull back the ear to stretch the auditory canal. When measuring ear temperature, insert the probe into the ear first (refer to the Picture 7.2), then press the button of scan. One second later you will hear "Beep", and you can remove the thermometer. The measurement is completed.

(Picture 7.1)

(Picture 7.2)

(Picture 7.3)

- 4) When the measured temperature exceeds 38.0°C the thermometer will send out "Beep-Beep-Beep" a long one and two short ones, and the display will show the icon of fever (Refer to the Picture 7.3).

Notes:

- Keep the probe surface clean, otherwise the measurement can't be accurate.
- Clean the probe with clean cloth or paper moistened with water or alcohol, and conduct measurement after the water on the surface of the probe all evaporates.

Taking Temperature in Forehead Temperature Mode


- 1) When the system is in the state of time / environmental temperature display, press the button of scan, and the last temperature reading will be displayed. Then, the system will go "Beep-Beep" and enter into the mode of measurement (Refer to the Picture 8.1-8.3).
- 2) Make sure the probe cap is on. Press the probe of the thermometer directly on the forehead.
- 3) Press the Scan button. One second later you will hear "Beep", and you can remove the thermometer. The measurement is completed.
- 4) When the measured temperature exceeds 38.0°C, the thermometer will send out "Beep-Beep-Beep", a long one and two short ones, and the display will show the icon of fever (Refer to the Picture 8.4).

(Picture 8.1)

(Picture 8.2)

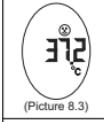(Picture 8.3)

(Picture 8.4)


80 English

Notes:


- When the product is used to measure forehead temperature, the reading is for reference only since skin temperature is significantly affected by the environment.
- When there is sweat or any other things on the forehead, e.g. cosmetics, the measured temperature can't be accurate.

Battery Power and Replacement

When the quantity of electric charge is not sufficient, the icon of low battery will appear on the lower part of the screen. Though the battery still can be used, it would be perfect to replace it (Refer to the Picture 9.1).

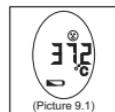
When the quantity of electric charge of the battery is lower than the minimal value, "Lo" symbols will be displayed on the screen, and battery icon flashes and the sound of "Beep" is sent out (Refer to the Picture 9.2). Unless you change a new battery, you can not use the thermometer to conduct measurement.

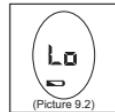Operation of changing battery

Take out the old battery: open the battery cap (Refer to the Picture 9.3). Put a new battery on: put on a CR2032 lithium battery with cathode downward into the battery seat (Refer to the Picture 9.4).

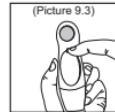Handling of Batteries

 **Warning:**

- Keep batteries away from children. Batteries are no toys.
- Batteries can be swallowed, which can be fatal. If a battery is swallowed, immediate medical assistance is crucial.
- Never throw batteries into fire. **Danger of explosion!**
- Do not force batteries open.
- If a battery leaks take care not to rub the battery acid into eyes or mucous membranes. Wash your hands after contact, flush your eyes with clear water and consult a physician in case symptoms persist.

Technical Specifications

Model	FT 4925 (FDIR-V1)
Measurement range body temperature/scan mode	32.0°C-42.9°C (89.6°F-109.2°F)
Measuring range room temperature	10.0 °C-40.0 °C (50.0 °F-104.0 °F)
Accuracy for body temperature measurement	±0.2°C/0.4 °F 35.5°C-42.0°C (95.9°F-107.6°F) ±0.3°C/0.5 °F 32.0°C-35.4°C (89.6°F -95.7°F) ±0.3°C/0.5 °F 42.1°C-42.9°C (107.8°F -109.2°F)
Resolution of display	0.1°C/0.1°F
Operating environment	10.0°C-40.0°C (50.0°F-104.0°F), Rh80% or below
Storage environment	-25.0°C-55.0°C (-13.0°F-131.0°F), Rh95% or below
Power supply	1x 3 V CR2032 lithium battery
Classification	II a (MD, Annex IX, rule 10)
Reference criterias and certifications according to EU directive for medical devices 93/42/EEC (MDD)	IEC/EN60601-1-2(EMC) CE0413 ISO 13485 (QM system)

82 English

For questions please contact:

Representative for Europe:

VISIOMED FRANCE
111 Rue De Stalingrad
93100 Montreuil, France
Phone: 0033-1-4808-44 13
e-mail: contact@visiomed-lab.com

Disposal of batteries

Your device is operated by batteries. When these batteries are exhausted, adhere to the following instructions:

△ Caution:

- Do not dispose of batteries with your normal household waste!
 - As user you are obligated by law to return exhausted batteries.
-
- The batteries must be removed from the device.
 - You may return your old batteries free of charge at the place of purchase.
 - Contact your local authorities for information on public collection points.

Environmental protection

Do not dispose of this device with the normal household waste.
Return it to public collection points for recycling. By doing so,
you will contribute to environmental protection.

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